# PHILIPS

**M1920A, M1921A, M1923A, 989803159571,  
989803159581, 989803160241, 989803160251,  
989803160261**

**Microstream® EtCO<sub>2</sub>**

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For more information, please call your local Philips sales office listed in your telephone directory or a Philips regional office listed below for the location of your nearest sales office.

## United States:

Philips Medical Systems  
North American Corp.  
3000 Minuteman Road  
Andover, MA 01810, USA  
(800) 225-0230

## Europe, Middle East and Africa:

Philips Medizin Systeme Böblingen GmbH  
Hewlett-Packard-Str. 2  
71034 Böblingen  
Germany  
Fax: (+49) 7031 463 1552

## Canada:

Philips Medical Systems Canada  
281 Hillmount Road  
Markham, ON  
L6C 2S3, Canada  
(800) 291-6743

## Asia Pacific Headquarters:

Philips Medical Systems  
30/F Hopewell Centre  
17 Kennedy Road  
Wanchai  
Hong Kong  
Tel: (852) 2821 5888  
Fax: (852) 2527 6727

## Latin America Headquarters:

Philips Medical Systems  
1550 Sawgrass Corporate Parkway #300  
Sunrise, FL 33323, USA  
Tel: (954) 835-2600  
Fax: (954) 835-2626

<http://shop.medical.philips.com>

Microstream®, FilterLine® and VitaLine™ are registered trademarks of Oridion Medical 1987 Ltd. ©2008 Oridion Medical 1987 Ltd. All rights reserved. This product or elements thereof is covered by one or more of the following US Patents: 6,422,240; 5,657,750; 5,857,461; 6,437,316; and 6,926,005 and their foreign equivalents. Additional patent applications pending.

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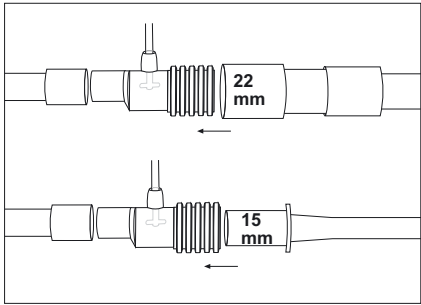
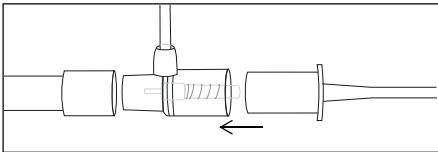
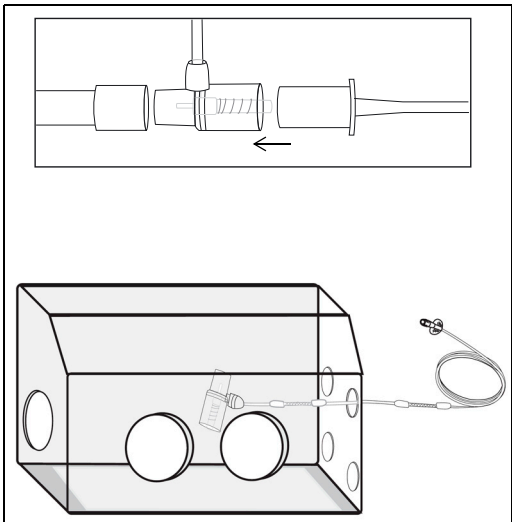
English

**Microstream® EtCO<sub>2</sub> Consumables for Intubated Patients**

These products are intended for use with the Philips M3015A Microstream® CO<sub>2</sub> Measurement Server Extension and the M8105A IntelliVue MP5 Monitor with Microstream CO<sub>2</sub>. On other Philips instruments they may be used only, if listed as accessories in the instrument’s Instructions for Use. Dispose of Microstream® products according to standard operating procedures or local regulations for the disposal of contaminated medical waste.

When using a VitaLine™ H Set Infant/Neonatal with an incubator, position one Nafion® humidity filter inside the incubator and the second Nafion® humidity filter outside the incubator, close to the incubator tubing access door.

Long FilterLines increase the sampling and alarm delays by 3 seconds compared to the specification for the normal FilterLines published in the instructions for use of the monitor.

**Ordering Information**

Product Brand Name	Patient	Product Numbers		
	Description	Standard (2m)	Longline (4m)	Vitaline
Filter-Line Set	Adult/ Pediatric	M1920A	989803160241	
FilterLine Set H	Adult/ Pediatric	M1921A	989803160251	989803159571
	Infant/ Neonatal	M1923A	989803160261	989803159581

 Do not flush the FilterLine® or VitaLine™ with any gas or liquid. This will damage the capnograph monitor and invalidate the warranty.

 These products do not contain latex.

**WARNING:** When used with a closed suction system, do not place the airway adapter between the suction catheter and the endotracheal tube. This is to ensure that the airway adapter does not interfere with the functioning of the suction catheter.

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